

type I, 5'-CTggCTgACTTATgCTTTTTACTgACT-3' (SEQ ID NO:1) and 5'-gATgCAggTgACTTTggCTACA-3' (SEQ ID NO:2)' (PCR product size 236 base pairs) and for angiotensin II receptor type II, 5'-ATTACTCCTTTTggCTACTCTTCCTC-3' (SEQ ID NO:3) and 5'-ggTCACGGGTTATCCTgTTCTTC-3' (SEQ ID NO:4) (PCR product size 489 base pairs). PCR amplifications were performed with 10 ng of template cDNA using a MJ Research PCR Cycle machine and the following PCR cycles: 1) 94 °C/2 min, 2) 94 °C/10 sec, 60 °C/30 sec, 72 °C/15 sec for 35 cycles, using High Fidelity Taq polymerase (Boehringer Mannheim) with components provided in the manufacturer's kit. Products of the PCR amplifications were identified by electrophoresis through a 0.8% agarose/TBE gel (Sambrook J, Fritsch EF, Maniatis T. Molecular Cloning: a Laboratory Manual, 2<sup>nd</sup> Ed., Cold Spring Harbor Laboratory Press, Cold Spring Harbor, NY, 1989). To confirm the identity of the PCR amplification products, the DNA was eluted from the gel and cloned in to the A/T cloning vector pMOSBlue (Amersham). Colonies containing a DNA insert of the correct size (Table 1) were fully sequenced on both strands to confirm their identity.--

---


#### REMARKS

Applicants have amended the Specification, in response to the Notice of Non-Compliant Amendment. This amendment is believed to be compliant with 37 CFR 1.1.21.

Attached hereto is a marked-up version of the changes made to the claims by the current amendment. The attached page is captioned **"Version with markings to show changes made."**

In view of the foregoing, Applicant submits the Application is now in condition for allowance and respectfully requests early notice to that effect.

Respectfully submitted,

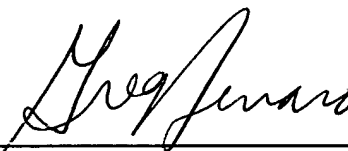

---

Gregory D. Ferraro  
Attorney for Applicants  
Reg. No. 36,134

Novartis Corporation  
Patent and Trademark Dept.  
564 Morris Avenue  
Summit, NJ 07901-1027  
(908) 522-6923

Date: May 6, 2002

Paragraph beginning at page 3, line 12 has been amended as follows:

PCR (Polymerase Chain Reaction) products were prepared and used as follows: Oligonucleotides specific to human angiotensin II receptor type I (GenBank accession number M93394) and angiotensin II receptor type II (GenBank accession number U15592) were designed using the Oligo 5.0 software programme to homologous regions from both sequences (Table 1). CDNA from human bone tissue was prepared following standard methods (Sambrook J, Fritsch EF, Maniatis T. Molecular Cloning: a Laboratory Manual, 2<sup>nd</sup> Ed., Cold Spring Harbor Laboratory Press, Cold Spring Harbor, NY 1989). For PCR the following oligonucleotide primer pairs were used: angiotensin II receptor type I, 5'-CTggCTgACTTATgCTTTTTACTgACT-3' (SEQ ID NO:1) and 5'-gATgCAggTgACTTTggCTACA-3' (SEQ ID NO:2) (PCR product size 236 base pairs) and for angiotensin II receptor type II, 5'-ATTTACTCCTTTTggCTACTCTTCCTC-3' (SEQ ID NO:3) and 5'-ggTCACGGGTTATCCTgTTCTTC-3' (SEQ ID NO:4) (PCR product size 489 base pairs). PCR amplifications were performed with 10 ng of template cDNA using a MJ Research PCR Cycle machine and the following PCR cycles: 1) 94 °C/2 min, 2) 94 °C/10 sec, 60 °C/30 sec, 72 °C/15 sec for 35 cycles, using High Fidelity Taq polymerase (Boehringer Mannheim) with components provided in the manufacturer's kit. Products of the PCR amplifications were identified by electrophoresis through a 0.8% agarose/TBE gel (Sambrook J, Fritsch EF, Maniatis T. Molecular Cloning: a Laboratory Manual, 2<sup>nd</sup> Ed., Cold Spring Harbor Laboratory Press, Cold Spring Harbor, NY, 1989). To confirm the identity of the PCR amplification products, the DNA was eluted from the gel and cloned in to the A/T cloning vector pMOSBlue (Amersham). Colonies containing a DNA insert of the correct size (Table 1) were fully sequenced on both strands to confirm their identity.